

# PROTEIN GEL STAINS


**Brilliant Blue R 250 Protein Staining Solution**

**Destain free G 250 Protein Staining Solution**

**Silver Staining Kit (low background)**

**Rapid Silver Staining Kit**

**Reversible Membrane Staining Solution**

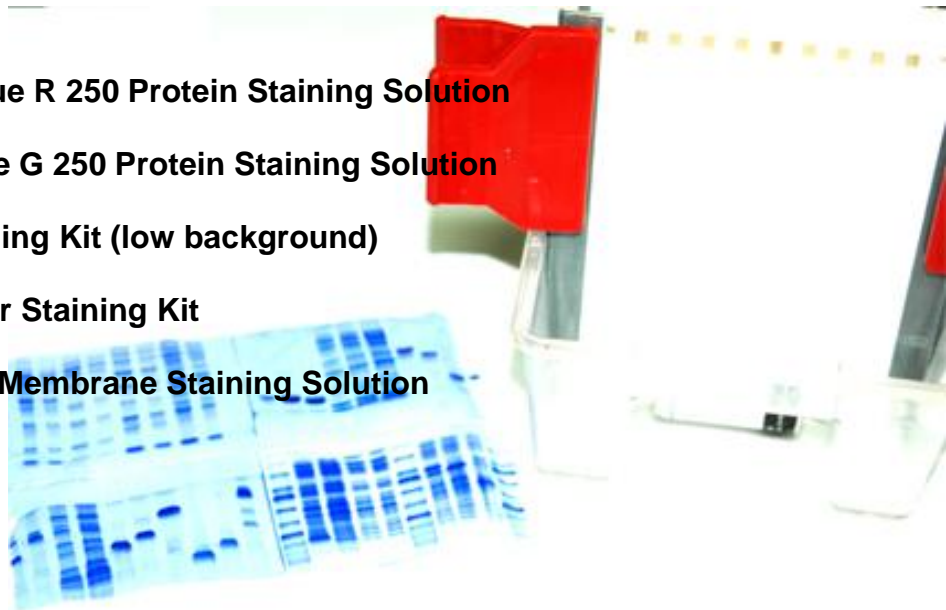

Staining Kits	Brilliant Blue R250 Staining Solution	Destain free G250 Staining Solution	Reversible Membrane Staining Solution	Rapid Silver Staining Kit	Silver Staining Kit (low background)
<b>Cat. No.</b>	EBP-1011	EBP-1021	EBP-1032	EBP-1052	EBP-1051
<b>Components</b>	1 L ready-to-use	1 L ready-to-use	1 L ready-to-use	500 ml SolA 500 ml SolB 500 ml SolC	500 ml SolA 500 ml SolB 500 ml SolC
<b>Protein Stain</b>	O	O	O	O	O
<b>DNA Stain</b>	X	X	X	O	O
<b>Sensitivity</b>	25-100 ng	5-10 ng with fixing 20 ng without fixing	100-200 ng	5-10 ng	1-5 ng
<b>Staining Time</b>	1 hr	> 2 hr	10 min	1 hr	>2 hr
<b>Background</b>	Low	Low	Low	Mid	Low
<b>Staining Step</b>	1. Staining (20min) 2. Destaining (30min)	1.Fixing (20min, optional) 2.Staining (1-4hr) 3.Destaining (DW, optional)	1.Staining (10min) 2.Destain (DW, optional)	1.Fixing (20min) 2.washing (10min) 3.Binding (20min) 4.Washing (1min) 5.Developing (10min)	1.Fixing (>20min) 2.Washing (20min) 3.Enhancing (1min) 4.washing (1min) 5.Binding (20min) 6.Washing (1min) 7.Developing (10min)